


# RHEOMETER REPORT

ABC

TEMP(C)	160.0
TIME(MIN)	4.0
TORQUE UNIT	kgf*cm

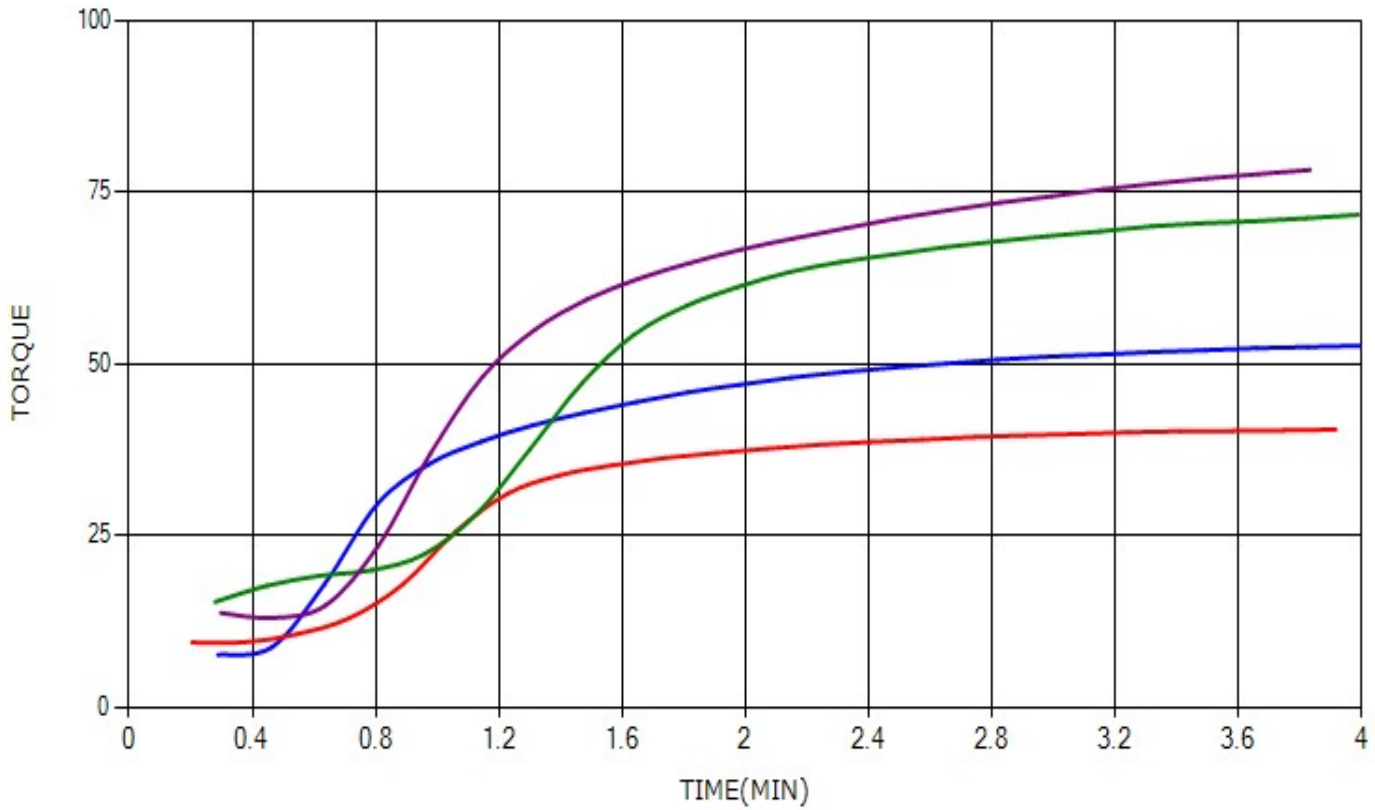

ID1	ID2	ID3	QC	ML	MH	ts1	ts2	tc10	tc50	tc90	Date
345			OK	8.8	52.7	0.63	0.63	0.63	0.97	2.32	2016/04/04 12:38:40
567			OK	9.5	40.5	0.54	0.71	0.71	1.06	2.07	2016/04/04 13:02:55
789			OK	17.8	71.7	0.62	0.79	1.13	1.47	2.65	2016/04/05 10:07:31
999			OK	13.1	78.3	0.64	0.81	0.81	1.15	2.66	2016/04/08 15:16:33

Lower Limit	8.0				0.50			
Upper Limit					1.50			
Average	12.3	60.8	0.61	0.74	0.82	1.16	2.42	
Standard deviation	4.09	17.36	0.04	0.08	0.22	0.22	0.28	
Coefficient of variation(%)	33.32	28.56	7.27	10.74	26.47	18.67	11.69	